
AKUBIONIC

ELECTROSTIMULATION
IONISATION

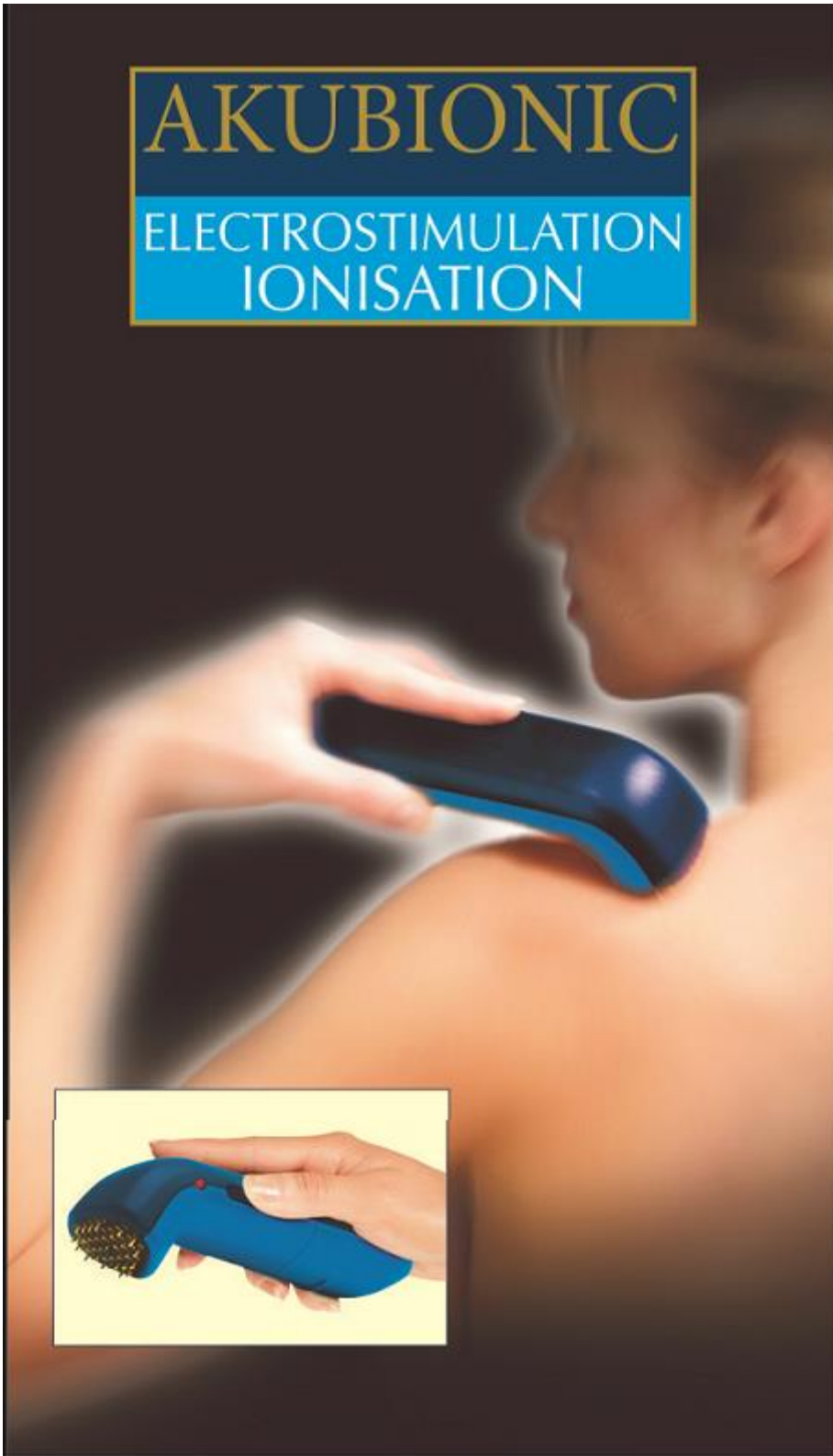AKUBIONIC

ELECTROSTIMULATION IONISATION

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MANUAL FOR USE

*HELPS IN EASING
ACHES AND PAINS AND
IN TONING UP*

AKUBIONIC

ELECTROSTIMULATION IONISATION

Akubionic's action reduces:

Contractures
Aches
Joint pain
Tendon pain

Akubionic's action increases:

Feeling of well-being
Vitality
Joint and muscle mobility
Remineralisation

Introduction

Akubionic was a great innovation in easing joint, muscle and tendon pains.

Now, Akubionic pioneers by combining electrostimulation and negative ions.

Akubionic is very practical and is used with only one hand ensuring you optimum well-being to your joints (knees, wrists, elbows, etc), lumbar and cervical regions, etc

Akubionic is highly compact and can go with you everywhere at all times.

It can be used at home, at work or on holiday.

How Akubionic works ?

Akubionic offers TWO exclusive actions to the body:

1) Electrostimulation of painful areas: with its 36 "electrostimulating" tips, **Akubionic** sends a flow of electrons to the sensitive area. These electrons are known for their soothing and easing effects.


2) Negative ion diffusion: there is a small carbon tip in the centre of the **Akubionic** head which sends a field of negative ions over the whole sensitive area.

A never-before feature!

• How to use AKUBIONIC ?

1) You simply switch on your Akubionic and apply to the area of the body that is troubling you. You will instantly feel slight prickling. The pulses are gradual and comfortable.

2) Please kindly note that the control function of the frequency and ionization was integrated to the intensity control in this latest version Akubionic.

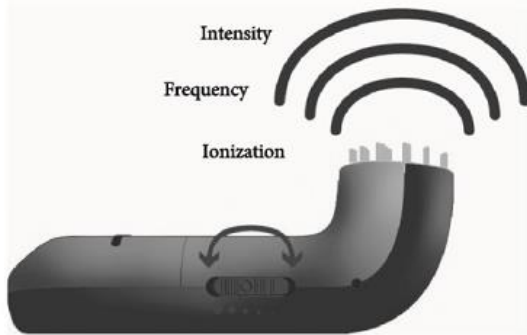The flows diffusion of electrons and negative ions is activated automatically when the device is switched on. The intensity modulation and the electrostimulation is done by pushing the unique wheel on the device.

Akubionic 2 is completely safety and natural

Precautions for use

Do not apply to device close to or onto injuries, wounds, lesions or rashes. Do not use if you have heart or coronary problems, if you suffer from thrombosis or in the event of pregnancy.

If in doubt, seek medical advice.

The following pages feature special detailed recommendations to ensure correct use of **Akubionic**.

Note: *this device must be kept out of children's reach.*

Using Akubionic :

Instructions for use:

To start the session, the following procedure must be adhered to:

Switch on the device. Turn the wheel on the side and first select the strength (start with the lowest).

Place Akubionic on the sensitive area of the body (avoiding the areas stated above) making small circular movements.

Then select the strength that suits you by turning the adjustment wheel on the device until you feel slight prickling.

This feeling should be pleasant. It indicates that you have reached the right strength.

You can change the strength during the session but it must remain tolerance.

A second button is provided to adjust the frequency of the electrostimulation flow and the negative ion diffusion according to how sensitive you are.

RECOMMENDATION :

For extra precaution, start off by selecting the lowest strength and frequency. When you use the highest strength and frequency, NEVER keep the head of the device in one and the same spot. Always stimulate the area chosen by making small circular movements over the skin. Sessions last from 2 to 7 minutes (or more according to the adjustments and depending on how sensitive you are). Do not press Akubionic too hard against skin to prevent any skin irritation.

COMMENT: slight reddening may appear during or after the session but this will quickly disappear.

Back pains

The most common pain is in the lower back. It can be caused by a sudden movement, bad posture, tight muscle, etc.

Stimulate the areas shown in the diagram opposite with **Akubionic**.

Shoulder pains

Origin: tight muscle, intensive exercise, stress, etc.

Make circular movements to stimulate the areas shown below with your **Akubionic**, adhering to the conditions stated in the "Precautions for use".

Elbow pain

Figure 1

Origin: tight muscle, intensive exercise, stress, etc.

Make circular movements over the areas shown below with your Akubionic.

Figure 2

Fig. 1, 2, 3 - Stimulation of the elbow area.

Figure 3

Knee pains

Figure 1

Origin: tight muscle, intensive exercise, stress, etc.


Make circular movement to stimulate the areas shown below with your Akubionic adhering to the conditions stated in the "Precautions for use".


Figure 2

Fig. 1, 2 - Stimulation of the knee area.

The following diagram shows the conventional acupuncture points which could help you to find optimum well-being

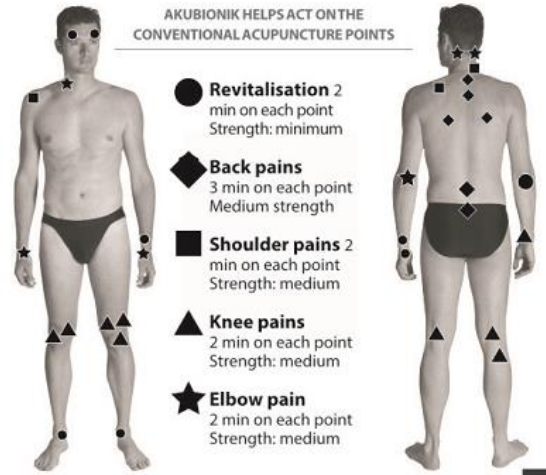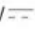1X6LR61 (not supplied) 9V 

WARNING: You should not dispose of this device with your household waste.

A selective collection system for this type of product is implemented by your local authorities. Please contact your local authorities to find out how and where collection takes place. These restrictions apply because electrical and electronic devices contain dangerous substances that have harmful effects on the environment or on human health and must be recycled.

This symbol indicates that electrical and electronic devices are collected selectively. The symbol shows a waste container crossed out with an X symbol.

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