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
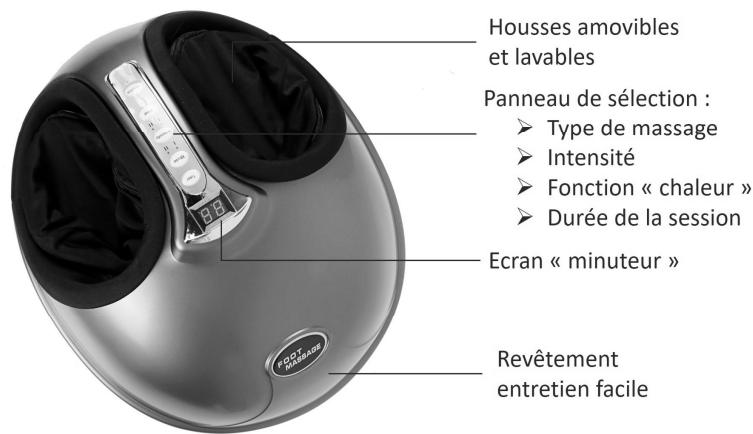
**Warning:** Please read this instruction manual carefully before using the instrument and keep it for future use.

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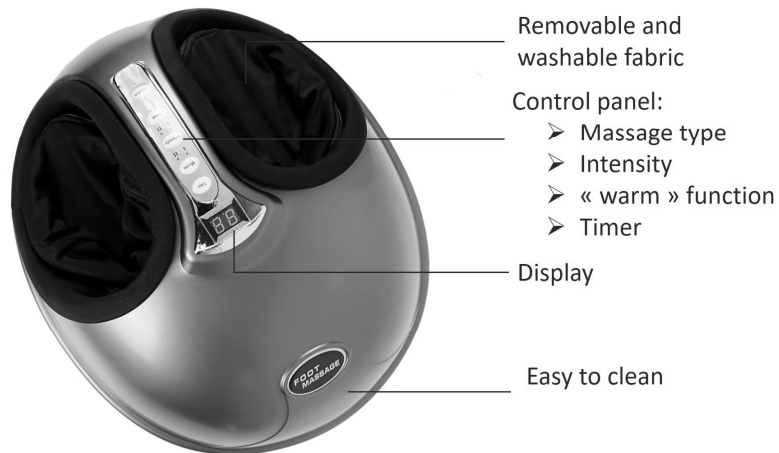Easy to clean

#### A. **GENERAL SAFETY INSTRUCTIONS**


1. The device is intended exclusively for private and not for commercial use. Only use this device as describe in this instruction manual.
2. Before to use, make sure that the supply voltage corresponds with the voltage marked on the rating label.
3. Do not place any part of this appliance in water or any other liquid.
4. Do not leave the device unattended when plugged in.
5. Do not handle the unit or plug with wet hands.
6. Unplug from the power outlet after use and allow the appliance to cool down completely before to clean it or to move it.
7. When disconnecting the cable, always pull at the plug; do not pull on the cable itself.
8. When using the device, always pull the power cord completely from its housing to prevent overheating.
9. If the cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid an electric shock.
10. Never try to repair it by yourself. It must be repaired by the manufacturer or its service agent or a similarly qualified person in order to avoid an electric shock.
11. The power cable must be checked regularly to avoid damage. If the power cable is damaged it must be replaced by a qualified service center to avoid any danger.
12. Always place your device on a dry surface.
13. Do not use outdoors.
14. Do not let the power cord on the edge of a table or work surface, and prevent it comes into contact with hot surfaces.
15. Keep it away from heat sources such as radiators, to avoid deformation of plastic parts.
16. Do not use coarse or abrasive sponges/cloths or steel wool to clean the appliance.
17. This unit is not designed to be operated by means of an external timer or a separate remote control system.

#### B. **SPECIFIC SAFETY INSTRUCTION**


1. Do not let water enter the product or be splashed near the surface of the product.
2. Do not operate the button with wet hands.
3. Keep the buttons from water or other liquid.
4. Water contact with electricity is dangerous. Do not leave the product in a humid environment (such as in the bathroom or shower room or swimming pool).
5. Before connecting the product to the power supply, check that the product voltage and supply voltage are the same.
6. If the plug, power cord or product itself is damaged, do not use the product.
7. In order to prevent the risk of electric shock, do not open the product outer cover.

8. Product is not intended for use by the physically or mentally impaired (including children), unless there is supervision or guidance by the people responsible for their safety when they use the product.
9. Children should not use this product. They should be under supervision to ensure that they will not treat this product as a toy.
10. If an unusual sound or noise is being heard from the product, turn off the product and immediately disconnect the power connection. Please contact Service center immediately.
11. Do not use this product in high temperature conditions or beside hot objects such as heaters, stoves.
12. Do not apply medical or cosmetic creams on the product or on your foot before using this product.
13. Do not destroy, bend stretch, twist or cut the power cord to avoid electric shock.
14. People who are sensitive to heat should not use the heat function of product.
15. *Recommended usage time for an individual is 15 minutes at a stretch.*
16. Please consult your doctor before using if you have the following situations
  - Receiving medical treatment or experiencing physical discomfort
  - Suffering from hypertension, heart disease, cerebral arteriosclerosis and cerebral thrombosis
  - Suffering from cancer
  - Pregnant or in a period
  - Suffering from osteoporosis
  - Using a pacemaker or implantable medical devices
  - Suffering from acute diseases
  - Minors
17. Do not fall asleep when using the product.
18. Do not use this product when you are drunk or feeling uncomfortable.
19. Do not directly contact the massager pad with your hands, feet joints, knee, chest and abdomen.
20. Do not use this product over thirty minutes in a row.
21. Do not use if the product's cloth cover has been damaged.
22. Do not use the product one hour before and after eating so as to avoid adverse effects.
23. Do not enable the heating function when the product is being used in a low room temperature. It is recommended that you should gradually adjust the temperature to abnormal range before enabling heating.
24. If the product is used immediately after being transferred from a cold environment to a warm environment, the accumulation of water vapor inside the product may affect the normal operation of the product or even cause mechanical failure. You are suggested to put the product in a normal room temperature for 1 hour before starting to use it after the transfer.

### C. PRODUCT DESCRIPTION
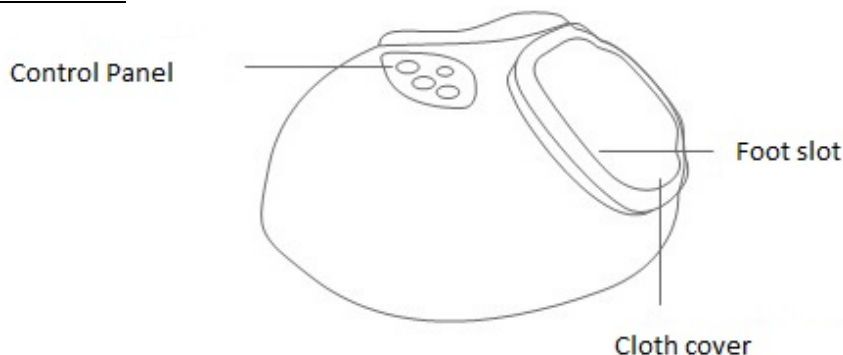
### Control Panel:

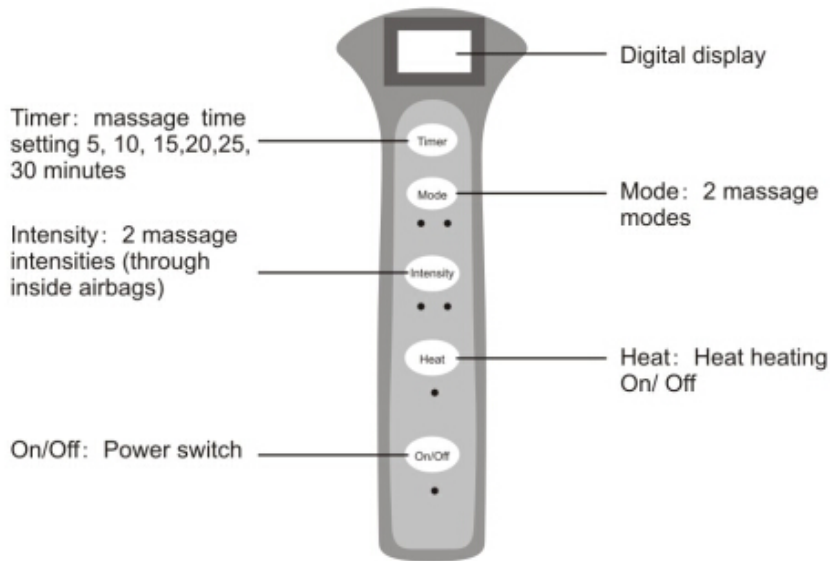- Untie the power cord and connect to power of 220V/50Hz.
- Press the Power Button ONCE, and the heating function is ON in a red indicator; the Intensity is low in a red indicator, and the product is on Mode 1 in a red indicator.
- Press the Power Button ONCE, then press the Heating Button ONCE, and the heating function will be OFF as the red indicator goes off.
- When the product is in operation, press the Intensity Button to switch between the two types of intensity. Red light indicates Intensity chosen.
- When the product is in operation, press the Mode Button to switch between two Modes. Red light indicates the Mode you have chosen.

### D. TROUBLE-SHOOTING

This page summarizes the most common problems you could encounter with the product :

Problem	Possible Cause	Suggestion
The product does not work	The product is not plugged in properly.	Ensure the voltage on the product corresponds to the local mains voltage and then plug in again.
	Power button not pushed	Push the power button
	The product has just session and has automatically switched to standby mode.	If you want another session, press the MODE button again.
	You may have used the product continuously for 30 minutes or more where the over heat protection has been activated.	Unplug the product and let it cool down for 60 minutes before you use it again.
	An excessive load has been applied on the product.	Reduce the excessive load and start again.

The product is abnormally warm	The product is likely to be overused.	Unplug the product and let it cool down before using it again.
Noise can be heard from the product	Noise is made by the mechanical parts of the product.	This is normal and does not mean that the product is faulty.

## E. CARE AND MAINTENANCE

### Cleaning

Do not use steel wool, abrasive cleaning agents or corrosive liquids (such as gasoline or acetone) to clean the product.

1. Unplug the plug from the power outlet.
2. Wipe with a dry cloth to clean the product.
3. Do not use a damp cloth or detergent
4. Do not hang the product with the cord.

### Storage

1. Unplug the plug from the power outlet.
2. Store the product in a cool dry place away from sunlight.
3. Do not put a heavy weight load on the product.

If the cover is damaged or there is any mechanical failure, please stop using the product.

## F. Specifications

Voltage 220V

Power 45W

*The European directive 2012/19/EC about electronic and electric waste, requires that you can't throw away defective domestic appliance with common waste. Used device has to be collect separately to optimize the recuperation rate and the recycling of the materials to reduce the impact on the health and the environment.*

*This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capacities, or by persons lacking in experience or knowledge, unless they are supervised or instructed on the use of this apparatus by a person responsible for their safety. Children should be supervised to ensure they do not play with this appliance.*

**Importé par / Importiert von / Geïmporteerd door / Imported by :**  
**CAPTELEC, 59170 Croix - FRANCE**

Carton et emballage papier à trier / Kartons und papierverpackung recyceln / Verpakkingsmateriaal in karton en papier te sorteren / Recycle carton and paper packaging