Gebrauchsanleitung | Instruction leaflet | Instruction d'utilisation
Instrucciones de empleo | Istruzioni per l'uso | Gebruiksaanwijzing
Руководство по эксплуатации | Instrukcja użytkowania | Bruksanvisning | Käyttöohje

promed KFS-100

Fußwärmer
Foot warmer
Chauffe-pieds
Scaldapiedi
Calentador de pies
Voetenwarmer
Грелка для ног
Ogrzewacz stóp
Fotvärmare
Jalkojenlämmitin

DE

EN

FR

IT

ES

NL


RU

PL

SE

FI

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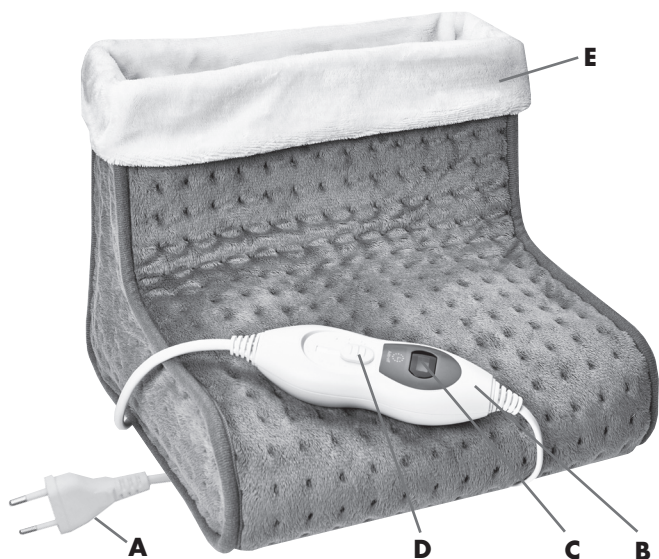- A. Netzstecker
- B. Bedienteil
- C. Display
- D. Temperaturregler
- E. Innenfutter

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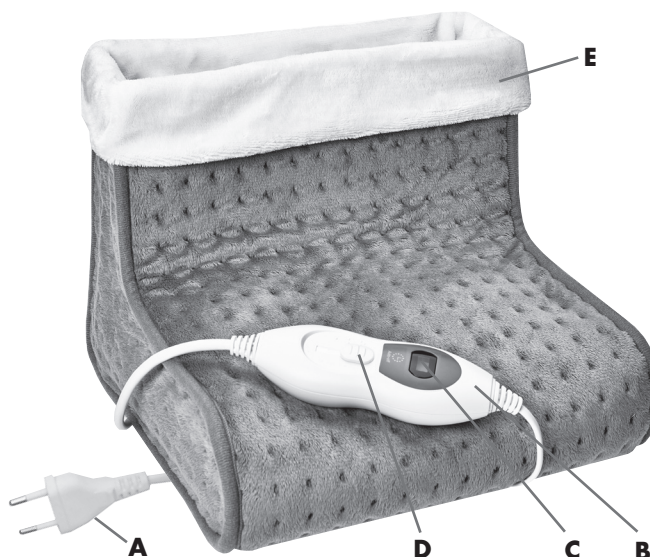
DESCRIPTION OF THE DEVICE

Please keep this user manual in a safe place!

Please read the instructions for use carefully before using the device for the first time.

Thank you very much for choosing the **Promed KFS-100** foot warmer. With this cosy, electric foot warmer, you will be able to comfortably warm up your feet. Using this device supports the blood circulation and thus the relaxation of tense muscle areas and increases general well-being after a stressful day.

The device is only intended for domestic use and not for commercial applications or use in hospitals or medical practices. To be able to make use of all advantages of your foot warmer, please read these instructions thoroughly and observe the care and safety instructions. This will allow you to enjoy the device for a long time. Please keep these instructions carefully for further use or if you wish to hand the device over to another person. We, the manufacturer, cannot be held liable for any personal injuries or damage to property resulting from the failure to comply with these instructions.

We hope you enjoy your new **Promed KFS-100**. With the information below, we would like to make you familiar with the device. Please read the instructions for use before using the device for the first time.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

SAFETY INSTRUCTIONS

Please note the following aspects prior to any commissioning:

- A.** The device may only be used with the controller supplied.
- B.** Voltage supply of the device: 220 - 240 V~.
- C.** If any problems occur with the device, please send it to the manufacturer.
- D.** Only for indoor use.
- E.** When using an electrical device near children, conscientious supervision is required.

DANGER!

- Never bring the device into direct contact with water or other liquids.
-   Never use the device in the bathtub, shower, pools, saunas or outdoors.
- Do not touch the device with wet hands.
- Only operate the device with the supply voltage specified.
- The device must not be used outdoors.
- This device uses supply voltage; therefore, the same safety regulations as for any other electrical device apply.
- Prior to any commissioning of the device, check that all components of the device are in perfect condition. In the event of a defect or a malfunction, immediately switch off the device and disconnect the power plug from power supply system.
- Any repair work may only be carried out by the manufacturer. Never open the device. The individual components of this device must not be disassembled, repaired or modified.
- Only use the device in dry rooms and if your body is dry.
- Never operate the device next to other sources of heat.
- Never carry or pull the device by holding it on the power cable. Disconnect the cable from the power socket using the plug and never using the cable.
- Do not pull the power cable over sharp edges or pinch it, do not let it hang down and do not place or suspended it over hot objects (e.g. hobs) or open flames and protect it against heat or oil.
- Only use an extension cable if you previously made sure that it is in a perfect condition.
- If the power cable of the device is damaged, have it immediately replaced by an expert.
- Never use the device after the power cable has been damaged, the device has been dropped or come into contact with water.
- Remove the power plug from the power socket if the device is not in use.
-  The device must not be used by infants and small children (0-3 years).
- This device can be used by children from the age of eight and people with limited physical, sensory or intellectual capabilities or with lack of experience and knowledge if they are supervised by a person responsible for their safety or have received extensive instructions in how to use the device by this person in such a manner that they understand the associated risks.
- Children must not play with the device.
- Cleaning and maintenance work must not be carried out by children without supervision.
-  An electrical device is not a toy. Children do not realise the risks which may result from using electrical devices. Therefore use and store the device out of reach of children.
- Disconnect the device from the power supply system before cleaning it. During cleaning, make sure that no liquids come into contact with the switch and that the device is only used in an absolutely dry condition.
-  Do not insert any needles or pointed objects into the device.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

SAFETY INSTRUCTIONS

WARNING!

- Never leave the device unattended if children or people unskilled in operating this device are around.
- Only use the device for applications described in these instructions for use.
- Never allow children to play with the packaging material, as there is a risk of suffocation.
- Do not place or use the device on hot surfaces or next to open gas flames to ensure that the housing does not melt.
- Do not store the device next to a washbasin or bathtub, since there is the risk that it can fall or be drawn into the washbasin or bathtub.
- Please consult a doctor before using the device if you suffer from medical conditions such as varicose veins, thrombosis etc.
- Never use the device on open wounds, burns, haematomata, swellings or similar injuries.
- Immediately stop using the device if it feels uncomfortable for you. Heat-sensitive people in particular should use the device carefully and only for short periods of application.
- The electric and magnetic fields generated by the device may interfere with the proper functioning of a cardiac pacemaker under certain circumstances. However, they do not exceed the limit values. Please consult your doctor and the manufacturer of your cardiac pacemaker before using the device.
Electric field strength: max. 5000 V/m,
Magnetic field strength: max. 80 A/m
Magnetic flux density: max. 0.1 mT
- Never use the device on animals.

DE

EN

Attention!

-  In order to ensure the proper functioning of the control element, do not place it on or underneath the device.
- The device must not be used by small children from the age of three unless the settings have been preset by a parent or legal guardian or unless the child has been sufficiently trained and instructed in how to operate the device.
- Do not use the device when it is wet.
- The device must not be used by helpless and defenceless people or people insensitive to heat.
- Children under the age of three must not use the device due to their inability to react to overheating.
- The electronic components in the switch may heat up when the device is used. Therefore, the switch must not be covered or lie on the device when it is in operation.

FR

IT

ES

NL

RU

COMMISSIONING THE DEVICE

PL

Before getting started, please check that the delivery is complete.

 Please use original parts only!

SE

 Under all circumstances, avoid accumulation of heat.

FI


COMMISSIONING THE DEVICE

1. To operate the device, insert the power plug **(A)** into the power socket.
2. In order to set the desired temperature, slide the temperature control **(D)** to the desired level (0 = OFF, 1 = low temperature, 2 = medium temperature 3 = maximum temperature). On the display screen **(C)**, the set level is now shown.
3. The device switches off automatically after approx. 90 minutes.
4. In order to restart the device after the operating time has expired, you must set the temperature control **(D)** to "0". After 5 seconds, you can start it by setting the temperature again.
5. After use, please disconnect the power plug from the power socket.

TROUBLESHOOTING

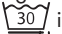Problem	Cause	Troubleshooting
The device does not warm up.	The device is not connected to the power supply system	Insert the power plug into the power socket and switch on the device
	The device is not switched on	Switch on the device
	The device is defective	Send the device to the manufacturer

CARE INSTRUCTIONS

- Before cleaning the device, the power plug must be disconnected from the power socket.
- Inner layer:  it can be washed by hand or in the washing machine at 30° C using a mild detergent.
- Foot warmer:  The device can be washed only by hand.
-  Do not iron the device.
-  Please do not use any detergents containing chlorine.
-  Do not dry clean the device.
-  The device and its cover cannot be tumble dried.
- Do not use any high-pressure cleaners.
- Make sure that the controller does not come into contact with any liquids. Liquids may result in damage.
- Never use the device while it is wet or damp. The device must not be switched on while it is drying.
- Let the device air dry. Do not use any auxiliary means (such as clothes pegs) to hang it up to dry. This may damage the cables inside.

Storage and regular inspection

- Do not crumple the device. Carefully fold the product to store it.
- Store the device in a dry place when it is not in use.
- Allow the device to cool before you fold or roll it up.
- Inspect the device and the power cable for possible signs of wear or damage at regular intervals.
- If there are such signs, if the appliance has been misused or does not work, return it to the supplier before switching it on again.
- Please do not use any chemical moth protection.
- Do not place any heavy objects on the device during storage, as this may damage the electrical wiring.
- In the event of malfunction or defects, send the device to the manufacturer.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

DISPOSAL

Power tools, accessories, and packaging should be recycled in an environmentally suitable manner.

Only for EU countries:

Do not dispose of power tools with the domestic waste! In accordance with the European Directive 2012/19/EU on waste electrical and electronic equipment and its implementation into national law, nonfunctioning power tools must be collected separately and recycled in an environmentally suitable manner.

This symbol indicates that disposal of this product in the domestic waste is not permitted within the EU. Waste devices contain valuable recyclable material that should be recycled. Recycling also prevents uncontrolled waste disposal from damaging the environment and human health. Therefore, please dispose of all waste devices in the appropriate collection systems or send the device to the place where you bought the device for disposal. They will then recycle the device.

WARRANTY

This product was produced with the utmost care and was thoroughly tested before it left the factory. Therefore, we provide a warranty period of 24 months from the date of purchase for this product. **Promed** products conform to their description and the relevant specifications; it is your responsibility to ensure that the product which you have purchased is designed for your intended use.

- In the event of material defects and/or manufacturing errors which appear during proper use and are detected during the warranty period, we will replace all defective parts of the product free of charge within the warranty period, including the labour costs for warranty repairs.
 - **The warranty does not cover:**
 - Normal wear on the product
 - Defects which were caused during transportation or storage of the product
 - Defects or damage which were caused through improper use or inadequate maintenance
 - Damage resulting from disregarding the notes in the operating manual
 - Damage resulting from modifications to the product which were not carried out by **Promed**.
 - Damage caused by sharp objects or as a result of torsion, compression, being dropped, an abnormal impact or other actions which are beyond **Promed**'s reasonable control.
 - Parts which are subject to wear (e.g. moving parts such as ball bearings, etc., fasteners, etc.) are generally excluded from the warranty.
 - **A claim under the warranty is not enforceable if:**
 - the product is not sent back in its original packaging or in adequately secure packaging;
 - it has been modified or repaired by a person or company other than **Promed** or one of **Promed**'s authorised dealers;
 - the product has been repaired using parts which are not approved by **Promed**;
 - the serial number / LOT number of the product has been removed, deleted, changed or made illegible.
 - For reasons of hygiene, any products which have been exposed to direct contact with a body or with bodily fluids (e.g. blood) must be packaged in an additional plastic bag before they are sent back. In this case, a specific note of this situation must be made in the package or in the accompanying documents.
 - For components or products which are replaced during the warranty period, a warranty is only provided for the remainder of the original warranty period; provided that the product or part replacement was carried out by **Promed** or one of **Promed**'s authorised dealers.
 - The warranty period begins on the day of purchase. Claims under the warranty must be established within the warranty period. Complaints arising after the expiry of the warranty period cannot be considered. Under these warranty conditions, the warranty shall only apply if the date of purchase can be proven by means of a sales receipt or similar.
- We reserve the right to make technical and visual changes as well as changes to the features.

DE

EN

FR

IT

ES

NL

RU

PL

SE

FI

WARRANTY

- This warranty is only valid and enforceable in the country in which it was purchased by the original buyer, provided that it was **Promed**'s intention that the product should be offered for sale in this country. This warranty is also enforceable in any country within the European Economic Area in which **Promed** has an authorised importer or sales partner. Depending on each country, specific and differing warranties and guarantees may be applicable owing to the particular applicable legislation. These legal provisions are neither excluded nor limited by these warranty conditions. To the extent permissible under national law, the warranty period shall not be extended, renewed or otherwise affected by any subsequent resale, repair or replacement of the product.
- The provisions of the UN Convention on Contracts for the International Sale of Goods shall not apply.
- Our warranty conditions do not affect the seller's legal warranty obligation.
- To the fullest extent possible under the applicable mandatory legislation, these warranty conditions constitute your sole and exclusive legal remedy and shall apply in place of any other express or implied warranty conditions. **Promed** is not liable for any unusual, incidental, punitive or consequential damages including, but not limited to loss of profits, loss of use, loss of revenue, costs of replacement equipment or facilities, third party insurance claims, and damage to property which are the result of the purchase or use of the product or occur owing to a breach of warranty, a breach of contract, negligence, product defects or any other legal or statutory circumstances, even though **Promed** was aware of the possibility of such damages. **Promed** is not liable for any delay in warranty claims.
- **Promed** cannot be held liable for any possible translation errors.

The following information is essential for smooth processing:

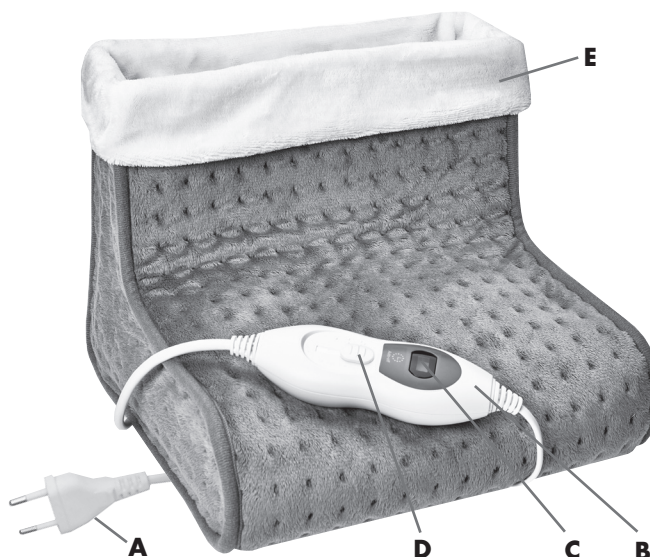
1. Original proof of purchase/receipt or dealer's stamp with date of purchase
2. Defect discovered
3. Product designation / type / serial/batch number

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Promed GmbH