EMS RELAXPAD

All languages available at: www.prorelax.com

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Manual de instrucciones

ES

GB Instructions for use

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IT

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NL

german
company
since 2009

Das CoolFit EMS RelaxPad

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
A Good Tool for Blood Circulation

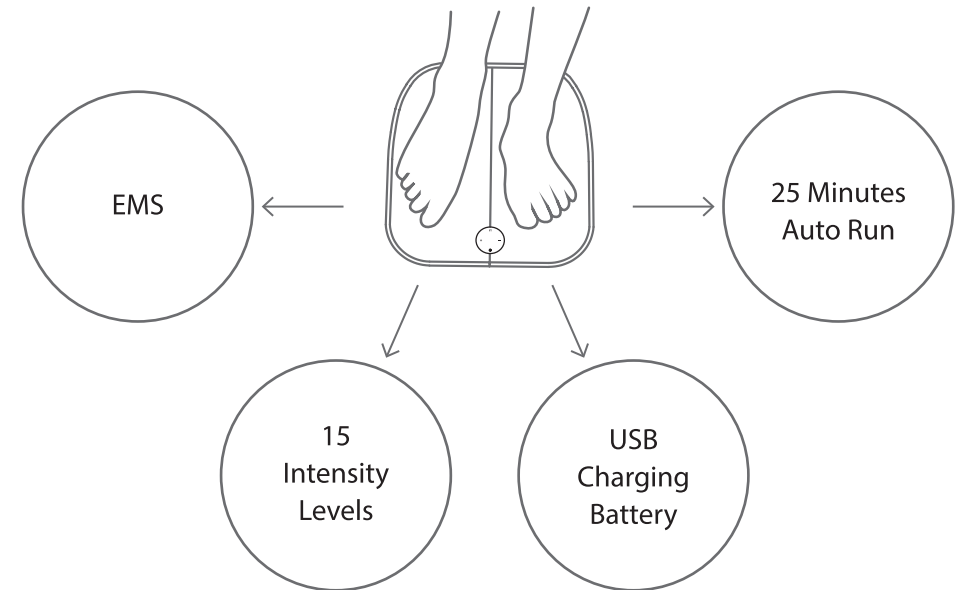Suggest one time per day.

Automatic turn off about 30 seconds later if DO NOT attach to Feet.

What's EMS Technology?

Electrical Muscle Stimulation (EMS), also known as neuromuscular electrical stimulation (NM ES) is the elicitation of muscle contraction using electric impulses. The impulses mimic the action potential coming from the central nervous system, causing the muscles to contract. EMS is proven to be quite effective as a proactive tool for muscle strengthening.

What is EMS used for?


- Muscle strengthening
- Rehabilitation purposes
- Preventing disuse muscle atrophy
- Increasing local blood circulation
- Relaxation of muscles
- Cosmetic muscle toning

Important

USE ONLY AS DIRECTED. THIS IS NOT A MEDICAL DEVICE. THE INFORMATION PROVIDED IS FOR INFORMATIONAL PURPOSES ONLY AND IS NOT INTENDED AS A SUBSTITUTE FOR ADVICE FROM YOUR PHYSICIAN OR HEALTHCARE PROVIDER

Know your CoolFit EMS RelaxPad


Main Unit

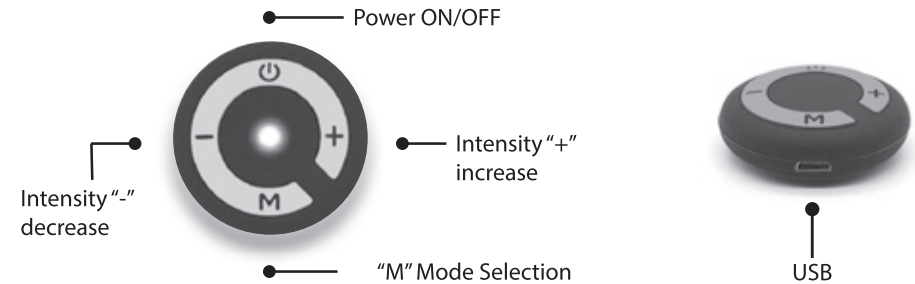Foot Mat


Remote control


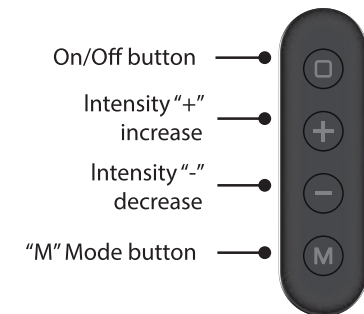
How to put a button cell battery CR2032 into the remote control

1. Pull out the back cover.
2. Put the battery inside.
3. Close the back cover.

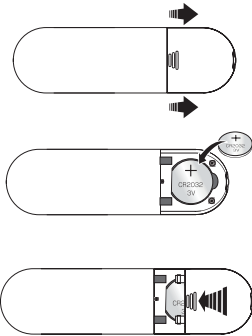Charge your EMS Foot Massager

- Plug in the main unit through micro USB charging cable (included)
- Connect the USB charging cable with AC/DC adaptor (not included)
- When charging completed, LED light will switch to white.

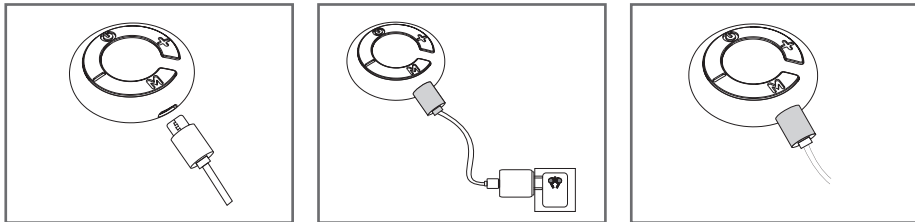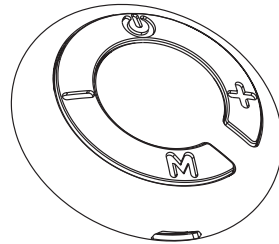How to operate

Preparations:

- Please don't wear socks or shoes when using the foot massager. Otherwise, you will not feel electrical pulses from your feet.
- Make sure both your bare feet must put separately on left and right side of the foot mat when power is on.
- It's better to wash or moisture your feet before using the foot massager. Dry feet may not feel electrical pulses due to high resistance.

Operation:


1. Put the main unit onto the foot mat, then make sure your bare feet are fully placed on the foot mat at each side.
2. Press the Power ON/OFF on the main unit to switch on the product.

Notes:

- When the product is powered on, the main unit starts from Level 1 in Mode 1.
 - When your bare feet take off from the footmat in working mode, the intensity level will be automatically down to Level 1.
3. Suggest one time per day.
 4. The remote control allows you to sit back and relax without having to bend over to operate the device.


5. In order to adjust intensity levels you can press the (+) button to increase or press the (-) button to decrease intensity. Intensity levels range from 1 to 15.

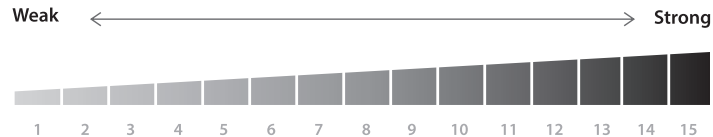6. Press M button to switch EMS Massaging Mode. LED light turns to green from white.

Auto Massaging Programme.

This 25 minute programme progressively increases massaging intensity in stage as time elapses. It automatically stops when the massaging is completed.

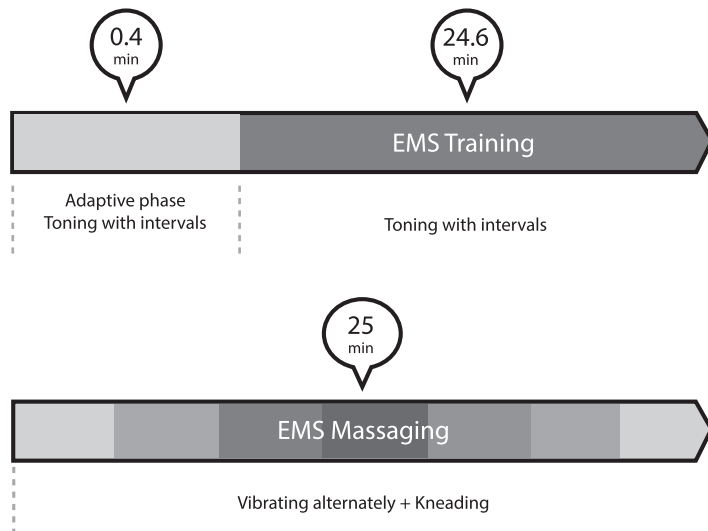General Safety Information

This EMS unit should never be used to mask or relieve undiagnosed pain. Before treating any symptoms of pain, you should consult your physician. For your safe and correct use of the device, it is recommended you to read the following guides before using it:

The following people must never use the foot massager.

- Anyone fitted with a pace maker or automatic implanted cardiac defibrillator (AICD)
- Anyone that is pregnant
- Anyone with deep vein thrombosis(DVT) to avoid mobilising the clot.

Precaution

Use a low intensity and shorter duration to introduce the muscles and avoid stimulation.

Other important safety precautions

1. Only use this device for the purpose of physical therapy.
2. Only use the accessories supplied with the foot massager.
3. Do not use this device in places with high humidity such as the bathroom.
4. Do not overload the electrical outlet.
5. If you are not going to be using the massager for a long time, remove the batteries from the remote controller. If the device malfunctions, disconnect it from the power source and contact customer services. Do not attempt to repair the device by yourself. As do so will invalidate the warranty.
6. Check the batteries in the remote control on a regular basis to ensure that they are operative.
7. Do not let the device come in contact with any metal object, such as belt buckle or jewelry.

LED-Anzeigestatus

Function	Operation	LED Indication Status
Power ON/OFF	Press 	ON/OFF
Increase +	Press 	Blink one time
Decrease -	Press 	Blink one time
Mode Selection	Press 	Color change (White/Green)
Charging	/	Red
Charging Completion	/	White
Non-Loaded	/	Blinking

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Function	EMS Foot Massaging
Power Supply	USB Rechargeable Battery, 3.7V
Frequency	4-35Hz/4-41Hz (with frequency change)
Intensity Level	15
Auto Shut-off	25 Minutes
Main Unit Size	Ø4.2cm x H 1.4cm
What's in the box?	1x Main Unit 1x Foot Mat 1x USB Charging Cable 1x Remote control 1x User Manual

WARRANTY PROVISIONS

Die KORONA electric GmbH, Willi-Bleicher-Str. 9, D-52353 Düren, (hereinafter called "KORONA") provides a warranty for the Prorelax® **RelaxPad** (hereinafter called the "Product") subject to the following preconditions and in the scope described below.

The warranty conditions described below do not affect the legally prescribed warranty obligations of the seller from the contract of purchase with the buyer. Moreover, the warranty applies irrespective of mandatory legal liability provisions, such as those according to the product liability law, or in cases of malafide intent or gross negligence, or owing to injury to life and limb, or to health by KORONA or its agents.

KORONA guarantees the defect-free operability and the completeness of the product, especially that the product is free of material defects, manufacturing defects and construction defects. The decisive factor in this regard is the state of the art of the science and technology at the time of manufacturing. The product must have manifested the defect that caused the damage at that time itself.

The warranty period is 2 years from the date of purchase of a new, unused product by the purchaser.

The warranty applies worldwide, i.e. regardless of the country in which KORONA has marketed the product and the purchaser has purchased it.

This warranty only applies to products that the purchaser has procured as a consumer and exclusively for personal purposes in the framework of domestic use (i.e. not for products that were procured or used in the framework of a commercial or independent professional activity).

If the product is found to be deficient in operability in accordance with the following provisions during the warranty period, KORONA will provide a free replacement according to these warranty conditions.

If the purchaser wishes to file a warranty claim, he must first contact KORONA customer service:

KORONA KUNDENSERVICE
Willi-Bleicher-Str. 9
Tel: +49 (0) 2933 90284 -80
D-52353 Düren
Germany
E-mail: service@korona-electric.de

The purchaser will then get further information on how to process the guarantee case, for example, where he can send the device free of cost and what documentation is required.

Filing a warranty claim can only come into consideration if the buyer can submit

- a copy of the invoice/proof of purchase and
- the original product

to KORONA or to a service centre nominated by KORONA upon demand.

Expressly excepted from this warranty are

- products that are used, cleaned, stored or maintained incorrectly and/or not in compliance with the conditions laid down in the operating manual, as well as products that were opened, repaired or modified by the purchaser or by a service centre not authorised by KORONA;
- products that were sold as second quality or used articles;
- consequential damages arising from a defect in the product (however, in this case, there may be entitlements from product liability or from other mandatory legal liability regulations).

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